XC-109

MyGuardian

QUICK START GUIDE

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Rev.1.0

GETTING STARTED

Thank you for purchasing XC-109 MyGuardian. Please follow the instructions and **read the safety instructions carefully to avoid any risk situations.**

ACTIVATE AND SET UP

1. Hold down Distress Button for **3 seconds**. Green LED light indicates the device is being activated.

2. Download and install "**XC-109 MyGuardian**" app to your Smartphone or you can scan QRcode below to finish the installation.

Package Contents

- XC-109 device
- Neck chain
- Battery included
- Quick start guide

System Requirements

- Android 4.4 or higher
- iOS 8 or higher
- Bluetooth 4.0

3. Follow **step-by-step** guide in the app to connect your mobile phone and XC-109 device.

The app will automatically scan device, please put XC-109 nearby your phone for following App settings.

The app will show the XC-109 device battery status after the device adds in. When the battery is lower than 15%, the app will automatically alert to the user.

Distress Contacts will automatically import from your mobile contact list and you can select distress contacts up to 6 after list imported.

IMPORTANT SAFETY INSTRUCTIONS

- **Wearing your XC-109 device at all times.**
- Always keep your smartphone around you.
- For correcting accuracy detection, please wearing your XC-109 device outside your shirt.
- Self send out distress SMS by **press Distress button at least 1.5 seconds.**
- Do Not water immersion the device more than 30 mins or up to depth 1 meter.

Typical of Bluetooth transmission distance is up to 10 meters, the device transmission range may vary depending on obstructions and your mobile performance, so please put your mobile phone close by.

Always wear XC-109 MyGuardian even in the shower or bath. The bathroom is often a place where people fall and get injured.

Beware of low battery status of XC-109 and your smartphone, it may be cause message transmission problem and getting risk of your life, therefore replace or refill the battery as soon as possible.

Highly Recommend **TOSHIBA CR2032** Lithium Battery, different brands may cause incorrect battery accuracy detection.

Warning! No fall detection technology is 100% accurate; if you fall and need help, don't wait for the XC-109 to activate, press the distress button immediately.

XC-109 MyGuardian with the Fall Detection pendant is intended to provide an extra level of protection. But please be aware that the Fall Detection pendant does not detect falls with 100% perfect accuracy or precision. The Fall Detection pendant is designed to detect falls based on the user falling down speed and impact force and therefore cannot detect every fall (for instance, if the user is gradually or slowly sliding down and/or without changing angle, especially from a sitting position, such as a wheel chair.), the user should always try to push the button if the user needs help.

DISCLAIMER AND NOTICE!

XC-109 MyGuardian will do whatever we can to support this installed device. We take responsibility for things in our control and will support you and the installation within our normal operating procedures. We cannot assume responsibility for the mobile network connection, cell phone towers or coverage, or other events or processes beyond our control.

XC-109 MyGuardian does not guarantee after send out distress SMS, the user can receive help or protection or that emergency help within a specific time frame.

TECHNICAL SUPPORT

For help with your XC-109 MyGuardian and mobile app installations, we provide illustrated instructions for set-up as well as resources for safety instructions. If our Quick Start Guide can not help your problem, please contact with our technical support at **service@honeywld.com.tw**.

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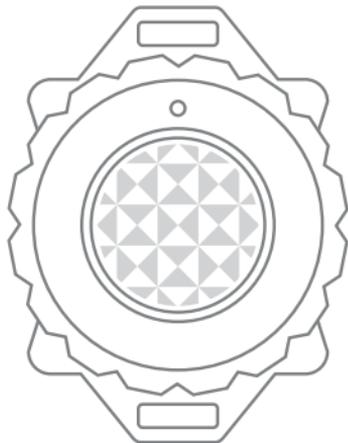
Model : XC-109 MyGuardian FCC ID : UNG-XC109

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Honeywld Technology Corp.

Made in Taiwan

Trademarks and Acknowledgements

All product and company names mentioned herein are the trademarks or registered trademarks of their respective owners. Any rights not expressly granted herein are reserved.

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Honeywld

Copyright © 2017 by Honeywld Technology Corp. All rights reserved