



Canturio® - Home Base Station



REF

43-5570-001-15

EDI



43557000115

QTY

1

SN

XXYYZZZZZZ

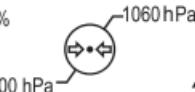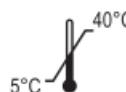

Rx
ONLY



IP22



FCC ID: 2AYAJ-BS1

CONTAINS FCC ID: 2AYAJ-WF1

CONTAINS FCC ID: R7T1201102



Use Only with Canary Medical
Supplied Power Source

MAC: A1:B2:C3:D4:E5:F6

Input: 5V d.c. === 0.1A Max

(01)00850015699153

(21)XXYYZZZZZZ



+H124435570001151B



+\$XXYYZZZZZD24BO

Distributed by:



Zimmer

1800 W. Center Street
Warsaw, Indiana 46580 USA

CANARY MEDICAL
2710 Loker Ave. West, Suite 350 | Carlsbad, CA 92010 USA
www.canarymedical.com | (833) 692-2627




| 8 | 7 | 6 | 5 | 4 | 3 | 2 | 1 |
|-----------|---|---|---|-----------|-----|-----------------|------------|
| | | | | REVISIONS | | | |
| | | | | ECO | REV | DESCRIPTION | AUTHOR |
| ECO101448 | | | | A0 | | Initial Release | C. Johnson |
| | | | | | | | 02/06/24 |

CONTROLLED DOCUMENT

CREATION TECHNOLOGIES


Feb 26, 2024

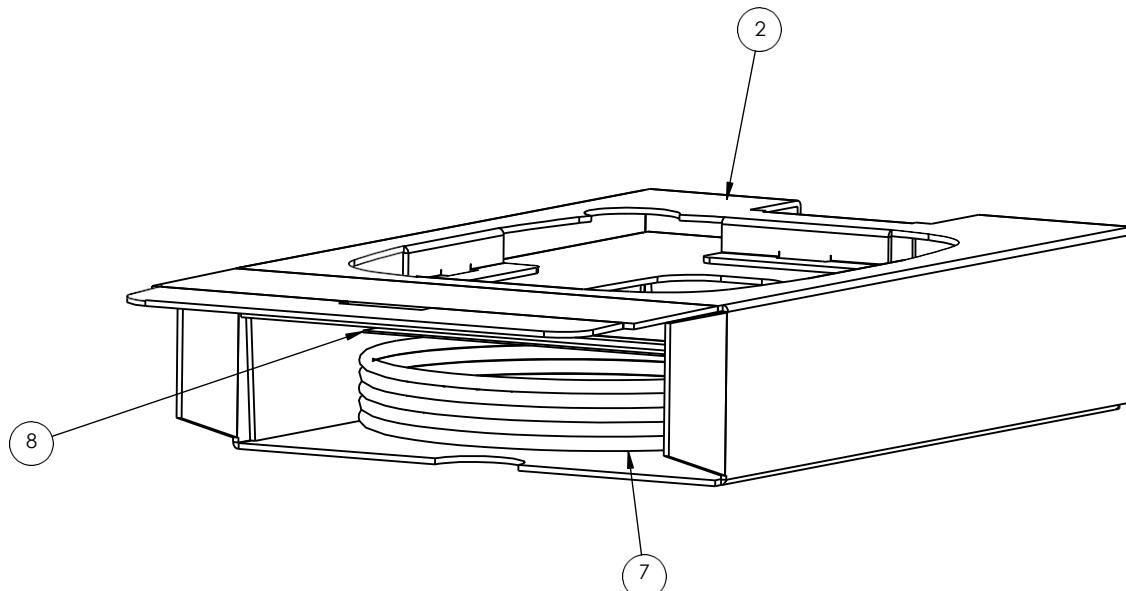

| ITEM NO. | Default/QTY. | PART NUMBER | DESCRIPTION |
|----------|--------------|-------------|---|
| 1 | 1 | 19005-0233 | TLA, Assy, Home Base Station, Gen 2 |
| 2 | 1 | 19005-0160 | Carton, Insert, Home |
| 3 | 1 | 19005-0210 | USB Power Supply |
| 4 | 1 | 19005-0161 | Carton, Outer, Bottom, Home |
| 5 | 1 | 19005-0139 | Carton, Outer, Top Home |
| 6 | 2 | 19005-0241 | Label, Shelf Pack, Home |
| 7 | 1 | 19005-0211 | USB Cable - 2m - Black |
| 8 | 1 | 19005-0243 | Quick Start Guide, Home |
| 9 | 1 | 19005-0159 | Bag, 10" x 8", clear, polyethylene, 2 mil thick, open-top |
| 10 | AR | 19005-0207 | Tape, Packaging |
| 11 | 1 | 19005-0224 | Box, Overshipper, Home Base Station |
| 12 | AR | 19005-0225 | Box, Multi-Pack, Home Base Station |

1. INSERT THE USB CABLE (ITEM 7) AND THE DOCUMENTATION (ITEMS 8) INTO THE LOWER OPENING OF THE INSERT (ITEM 2) AND FOLD TABS CLOSED.

| | | | |
|--|----------|------------|-------------------------------------|
| THIRD ANGLE PROJECTION | PROJECT | 19005 | Creation Technologies Design Center |
| | ENGINEER | C. JOHNSON | |
| | DATE | 2/5/2024 | |
| THE INFORMATION CONTAINED IN THIS DRAWING IS CONFIDENTIAL AND PROPRIETARY. ANY COPYING OR REPRODUCTION IN WHOLE OR PART WITHOUT WRITTEN CONSENT IS PROHIBITED. | | | |
| DESCRIPTION | | | |
| Assy, Packaged, Home Base Station, Gen 2 | | | |
| DRAWING NUMBER | SCALE | REV | SHEET |
| 19005-0236 | 1:25 | A0 | 1 of 4 |

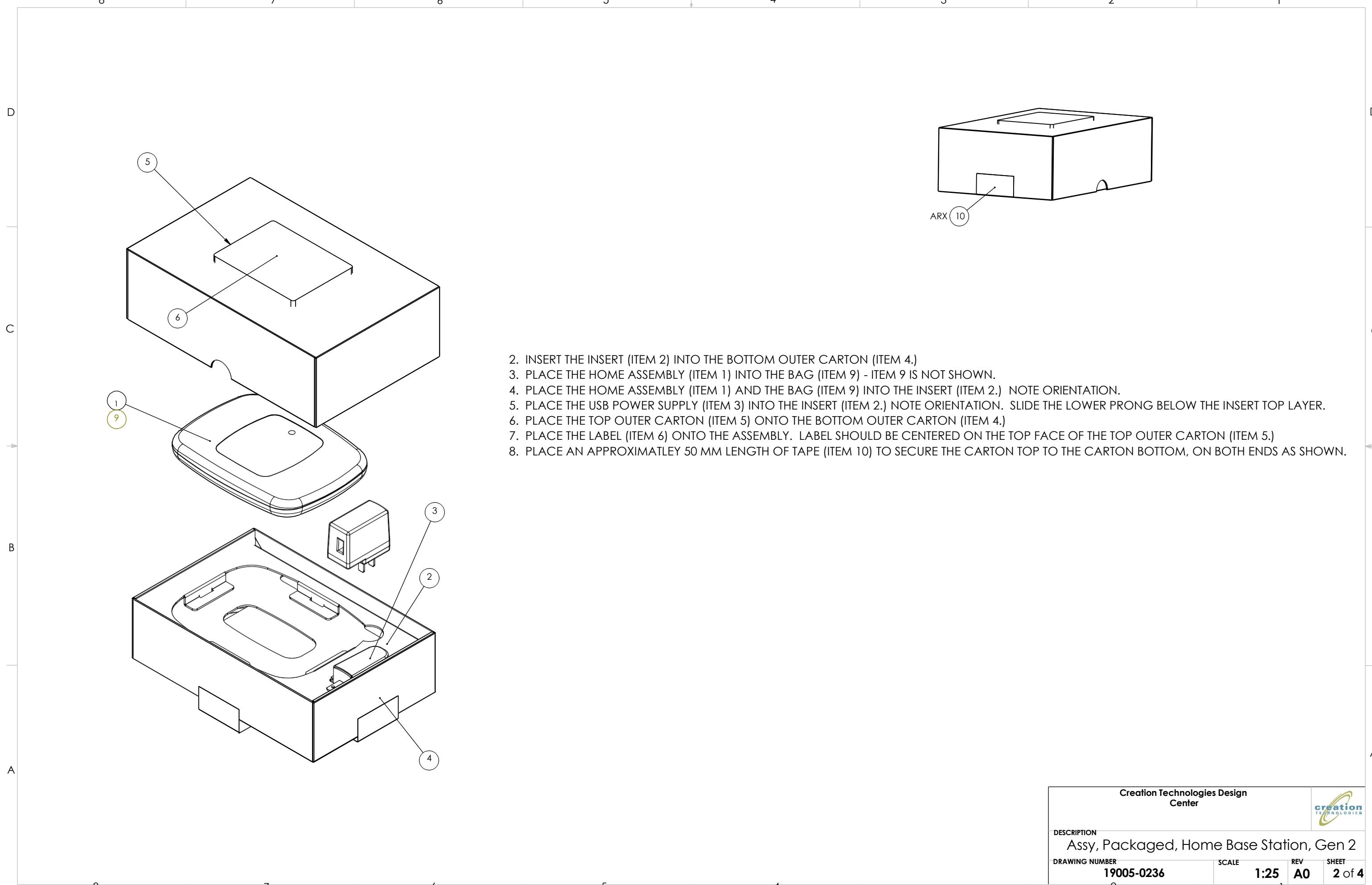

| | | |
|-------------------------------------|--|--------------------------|
| Creation Technologies Design Center | | creation TECHNOLOGIES |
| DESCRIPTION | Assy, Packaged, Home Base Station, Gen 2 | |
| DRAWING NUMBER | 19005-0236 | SCALE |
| | 1:25 | REV |
| | A0 | 2 of 4 |

8 7 6 5 4 3 2 1

D

D

C



C

B

B

A

A

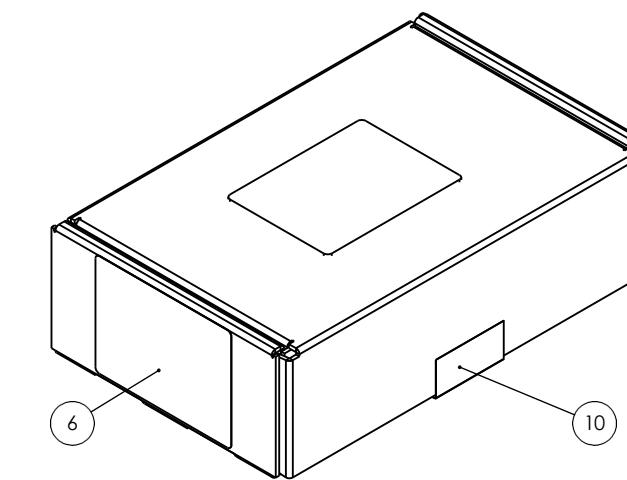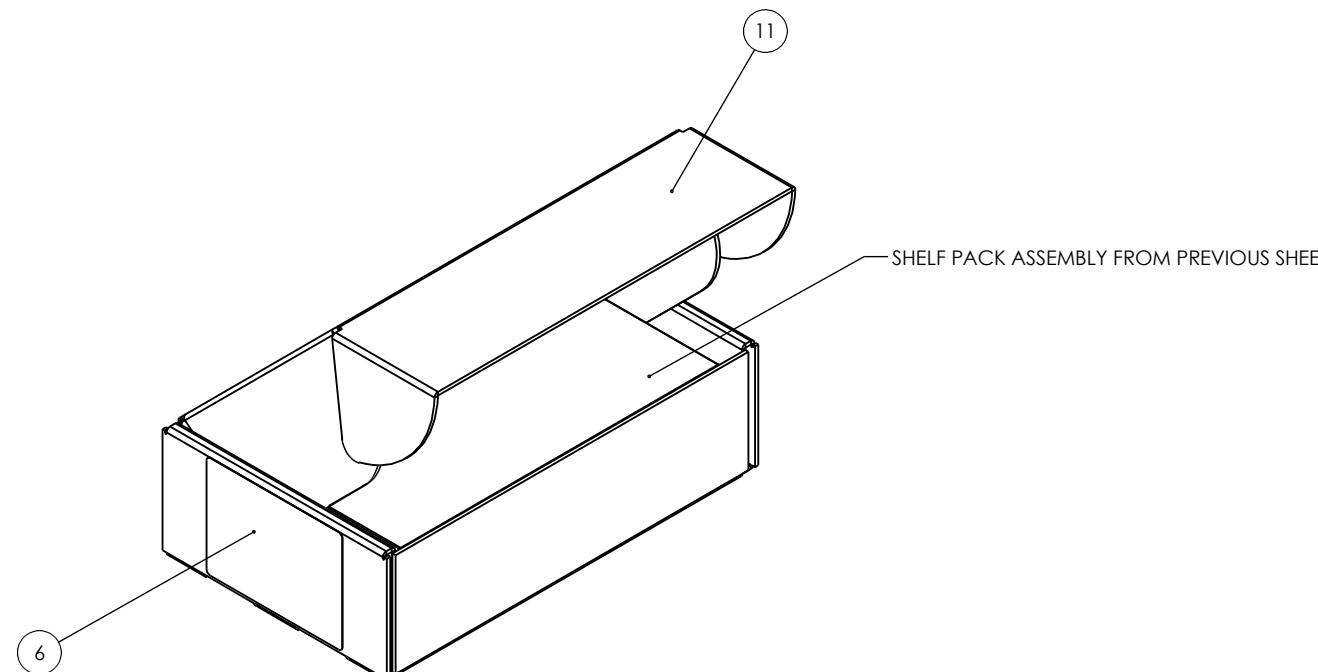

OVER-SHIPPER ASSEMBLY INSTRUCTIONS

9. PLACE THE SHELF PACK, FROM THE PREVIOUS SHEET, INTO THE OVER-SHIPPER (ITEM 11).
10. TAPE THE OVER-SHIPPER BOX SHUT USING PACKAGING TAPE (ITEM 10).
11. PLACE THE SECOND PACKAGING LABEL (ITEM 6) CENTERED ON THE SHORT SIDE OF THE OVER-SHIPPER BOX AS SHOWN.

8

7

6

5

4

3

2

1

| | | | |
|-------------------------------------|------------|---|--------|
| Creation Technologies Design Center | |  | |
| DESCRIPTION | | | |
| DRAWING NUMBER | 19005-0236 | SCALE | 1:33.3 |
| REV | A0 | SHEET | 3 of 4 |

8 7 6 5 4 3 2 1

D

D

C


C

B

B

A

A

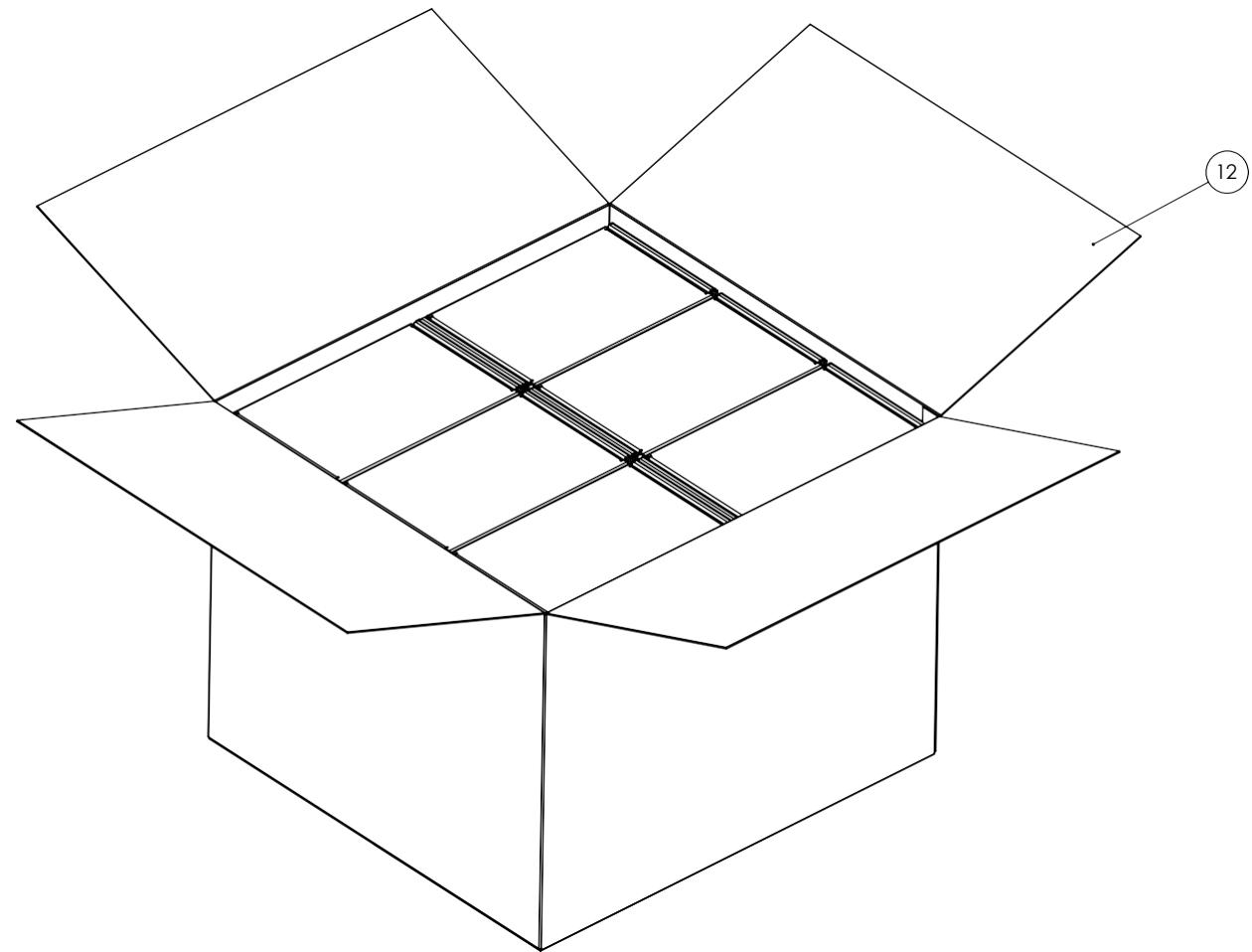

MULTI-PACK ASSEMBLY INSTRUCTIONS (OPTIONAL)

12. ERECT THE MULTI-PACK BOX (ITEM12) USING THE PACKING TAPE (ITEM 10) WHILE LEAVING THE TOP PORTION OF THE BOX OPEN.
13. PLACE THIRTY ASSEMBLED OVER-SHIPPER UNITS, FROM THE PREVIOUS PAGE, IN THE MULTI-PACK BOX.
14. SEAL THE MULTI-PACK BOX WITH PACKING TAPE.

| | | |
|-------------------------------------|-------|---|
| Creation Technologies Design Center | |  |
| DESCRIPTION | | |
| DRAWING NUMBER | SCALE | REV SHEET |
| 19005-0236 | 1:2 | A0 4 of 4 |

