



A Division of Instrumentarium

Made in USA, Manufactured by  
Datex-Ohmeda, Inc. Issaquah, WA


ULTRAVIEW<sup>®</sup> DIGITAL TELEMETRY  
ECG TRANSMITTER  
90341-05

FCC ID: CM676A90341-WMTS  
CANADA:

Operation requires prior coordination  
with the WMTS frequency coordinator

334-5553-01



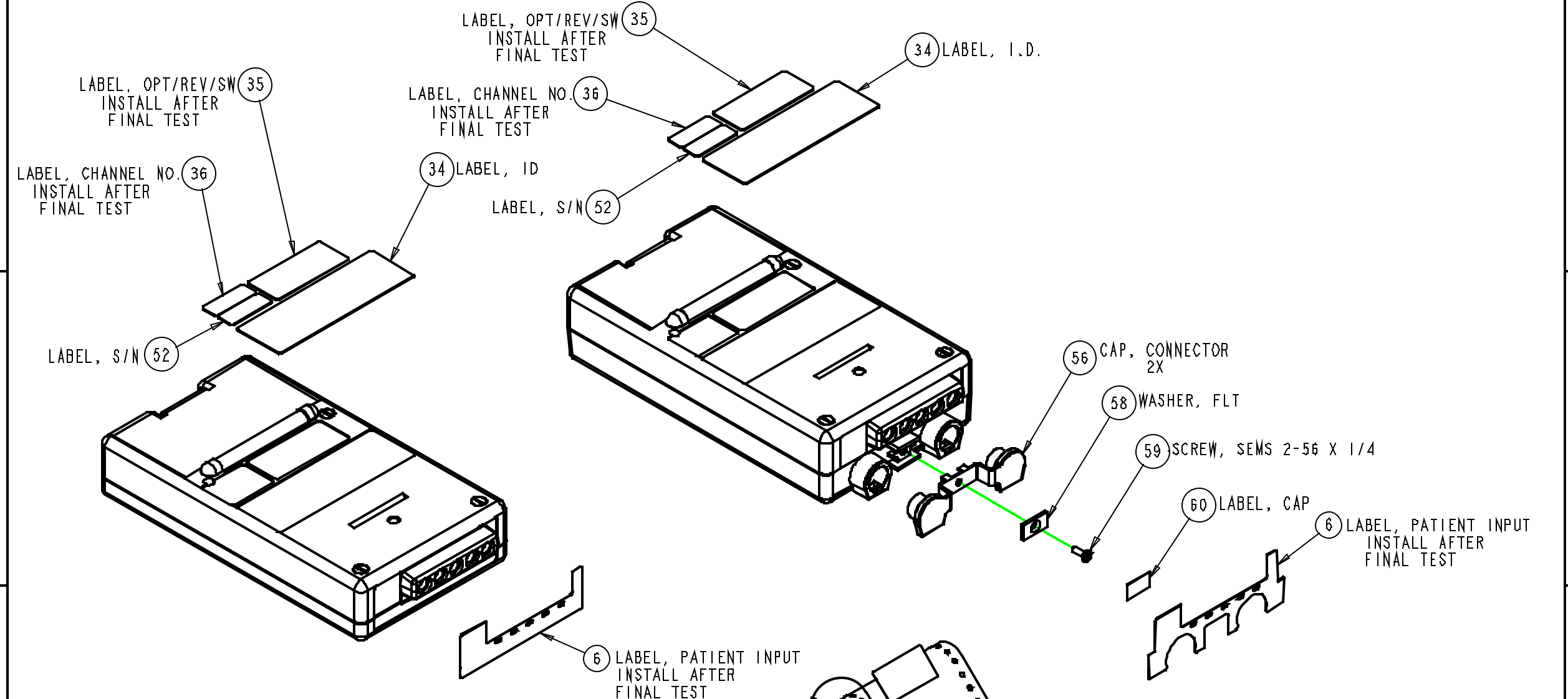

**ADHESIVE PREPARATION:**  
(MASTERBOND EP42HT)

ADHESIVE MUST BE MIXED BY WEIGHT  
INSTRUCTIONS FOR PREPARATION ARE AS FOLLOWS:  
PLACE CLEAN EMPTY ALUMINUM CUP TO BE USED FOR ADHESIVE ON  
GRAM SCALE. TARE THE SCALE, SUCH THAT DISPLAY READS 0.00  
WITH CUP ON SCALE.

MEASURE 10 GRAMS MINIMUM OF PART "A" USING THE GRAM SCALE,  
INTO THE ALUMINUM CUP.

MULTIPLY THE AMOUNT OF PART "A" X 0.4. THIS IS  
THE WEIGHT OF PART "B" TO BE USED WITH PART "A".

ADD THE NUMERIC VALUES OF ACTUAL WEIGHT PART "A" AND CALCULATED  
WEIGHT PART "B".

ADD PART "B" TO PART "A", IN THE ALUMINUM CUP ON THE SCALE,  
UNTIL THE TOTAL WEIGHT EQUALS THE AMOUNT CALCULATED AS  
"A"+"B" IN THE PREVIOUS STEP.

GELATION OF ADHESIVE IS APPROXIMATELY TWO HOURS.

**LABELING INSTALLATION**

CLEAN LABEL SURFACES WITH ISOPROPYL  
ALCOHOL. REMOVE LABEL FROM LINER  
HANDLE LABELS BY THE EDGES ONLY  
DO NOT TOUCH LABEL SURFACES OF CASE AND COVER.

IMPRINT LABELS ITEMS 35 AND 36 (WITH PERTAINING  
INFORMATION) USING THESE APPROVED TYPE WRITER  
RIBBON SOURCES:

O.S. EATON CORP. 99 MADISON AVE. NEW YORK, NY. 10016 (800) 821-2288 RIBBON: #172-B	MID CITY COLUMBIA 35 S ST. CLAIR ST. DAYTON, OH. 45401 (800) 462-2336 RIBBON: #CGL-79
--	---

SpaceLabs Medical, Inc. REDMOND, WASHINGTON		SIZE C	CODE ID NO 59873	PART NO 653-0095-00	REV L
DRAWN M. AMAYA 6-8-98	ISSUED	SCALE 1:1	FILE NO.	SHEET 10 OF 10	