


# ULTRAMATIC

E1004204  
REV E

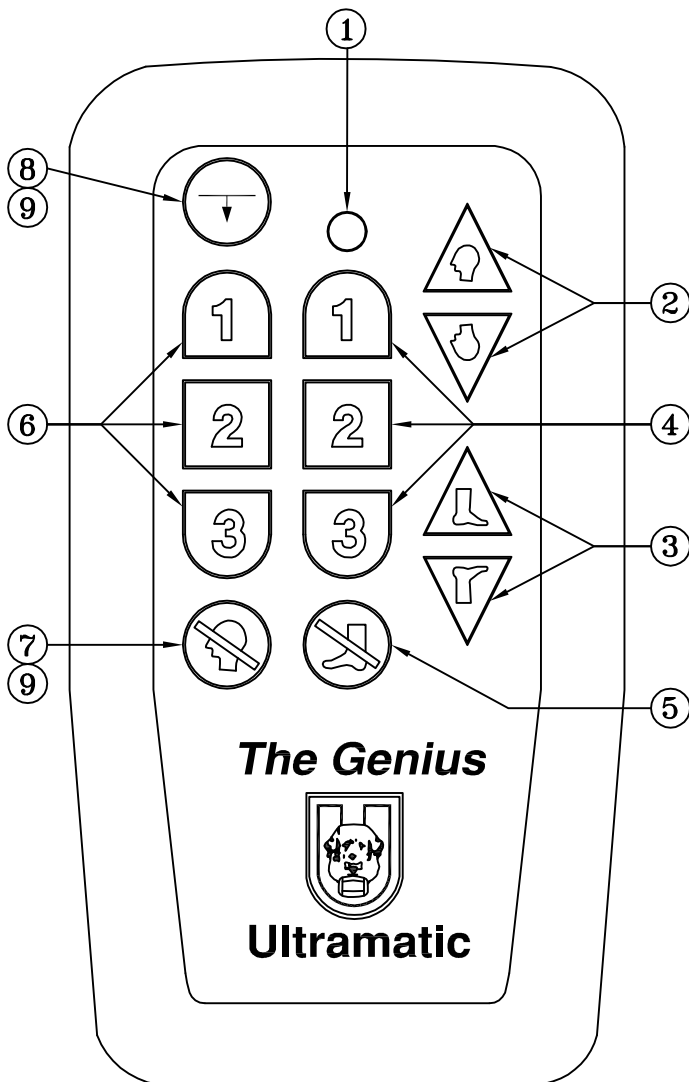

① **TRANSMISSION INDICATOR LIGHT:**  
Verifies a button is being pressed.

② **HEAD UP/DOWN BUTTONS:**  
Adjusts the head section of the bed to the desired position.

③ **FOOT UP/DOWN BUTTONS:**  
Adjusts the foot section of the bed to the desired position.

④ **FOOT MESSAGE ON/SPEED CONTROL BUTTONS:**  
Adjust message speed: 1 for LOW, 2 for MEDIUM, 3 for HIGH. After turned on, the message motor will shut off in 30 minutes.

⑤ **FOOT MESSAGE OFF BUTTON:**  
Turns off foot massage.

⑥ **HEAD MESSAGE ON/SPEED CONTROL BUTTONS:**  
Adjust message speed: 1 for LOW, 2 for MEDIUM, 3 for HIGH. After turned on, the message motor will shut off in 30 minutes.

⑦ **HEAD MESSAGE OFF BUTTON:**  
Turns off head massage.

⑧ **FLAT BUTTON:**  
Levels the head and foot sections to the flat position and turns off massage functions when button is pressed and released.


⑨ **AUTO LEARN FUNCTION:**  
Auto Learn is only active for 90 seconds once enabled. To enable Auto Learn, choose the option that meets the bed setup.  
• For Powerdown box equipped beds, press the Learn button on the Powerdown Box for 2 (two) seconds and release.  
• For beds without Powerdown, Unplug the bed from wall outlet for 20 seconds, then plug bed back into outlet to reset system.  
Next, simultaneously press and hold the flat and head massage off buttons. While holding these buttons the LED transmission indicator light (at the top of the wand) will turn off and then it will begin flashing rapidly. Continue to hold for 2 seconds then release. The new wand code has been learned. Please note that when using the auto learn feature when more than one wireless equipped bed is in transmission range (approx. 50 ft. radius) you must remove power from all beds not intended to learn the new hand control address.

---

## OTHER FEATURES:

- \* Both head and foot functions may be operated simultaneously.
- \* Radio Frequency is 310 MHz.
- \* To conserve battery life, the transmitter turns off all functions if any button is pressed for more than 50 consecutive seconds.
- \* **WARNING:** The transmitter and the receiver of the wireless control system have been tested and comply with FCC and Industry Canada rules. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

# WIRELESS RECEIVER

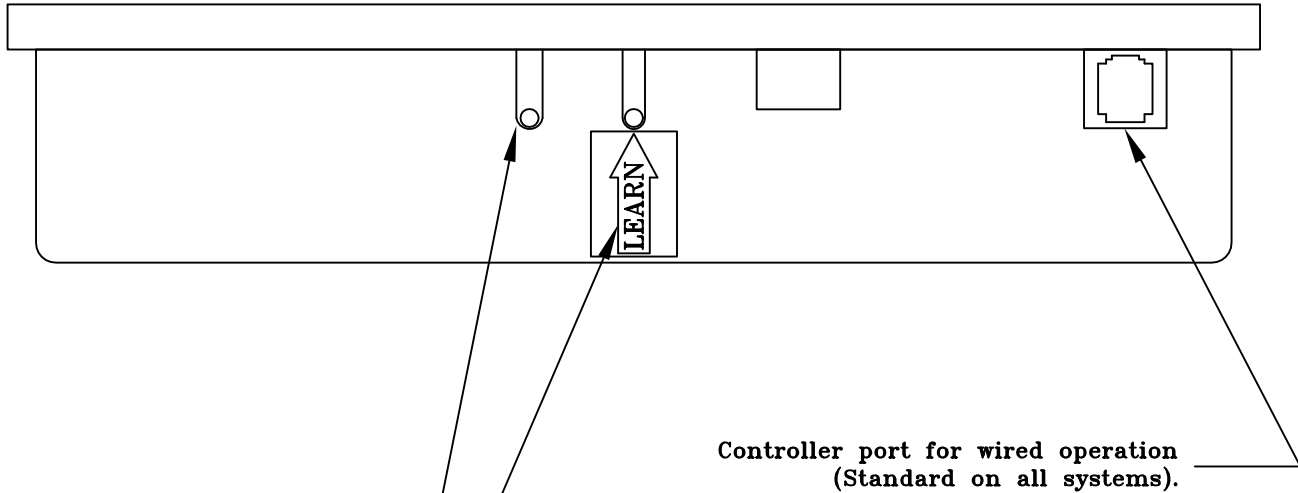

Controller port for wired operation  
(Standard on all systems).

## WIRELESS HAND WAND LEARN SWITCH:

To learn a new hand wand into the receiver, press this learn switch. Notice the learn light turns on and stays on until the receiver has learned its new hand wand or the learn period has ended (22 sec.). During the learn period press any button on the new hand wand. When the bed responds, the new hand wand is learned.

## LEARN/RECEIVE LIGHT:

This light turns on when the learn switch has been pressed and stays on to indicate the receiver is waiting for a hand wand to learn. When the receiver has learned a hand wand, this light will flash while the receiver is detecting a signal from any hand wand, learned or not.

---

---

## OTHER FEATURES:

- \* Each hand wand has its own unique address, programmed at the factory.
- \* Once a hand wand's address has been learned by the receiver, the receiver will respond to only that, most recently learned, hand wand.
- \* If power is lost (bed becomes unplugged, power outage, etc.), the receiver will remember its learned address, and pre-set bed positions (for systems with the programmable position feature).