



Electromagnetic Emission Compliance Report for FCC Certification

Test Report No. : E132R-002 AGR No. : A131A-174

Applicant : SD Biosensor, Inc.

Address : C-4th&5th Floor Digital Empire Building 980-3, Yeongtong-dong, Yeongtong-gu,

Suwon-si, Kyonggi-do, Korea

Manufacturer : SD Biosensor, Inc.

Address : C-4th&5th Floor Digital Empire Building 980-3, Yeongtong-dong, Yeongtong-gu,

Suwon-si, Kyonggi-do, Korea

Type of Equipment : SD GlucoNavii® Voice N Blood Glucose Meter

(Part 15 Class B Computing Device Peripheral)

Model Name : 01GM51

FCC ID. : RPJ01GM51

Serial number : N/A

Total page of Report : 11 pages (including this page)

Date of Incoming : January 23, 2013

Date of Issuing : February 04, 2013

SUMMARY

The equipment complies with the requirement of FCC CFR 47 PART 15 SUBPART B, Section 15.101.

This test report contains only the results of a single test of the sample supplied for the examination. It is not a general valid assessment of the features of the respective products of the mass-production.

Prepared by:

Seung-Hyun, Nam / Asst. Chief Engineer ONETECH Corp.


Approved by:

Gea-Won, Lee / Managing Director ONETECH Corp.

Report No.: E132R-002

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)

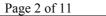



CONTENTS

	Page
1. VERIFICATION OF COMPLIANCE	4
2. GENERAL INFORMATION	5
2.1 PRODUCT DESCRIPTION	5
2.2 MODEL DIFFERENCES	5
2.3 RELATED SUBMITTAL(S) / GRANT(S)	5
2.4 TEST SYSTEM DETAILS	6
2.5 CABLE DESCRIPTION FOR THE TEST SYSTEM	6
2.6 TEST METHODOLOGY	6
2.7 TEST FACILITY	6
3. SYSTEM TEST CONFIGURATION	7
3.1 MODE OF OPERATION DURING THE TEST	7
3.2 EQUIPMENT MODIFICATIONS	7
3.3 CONFIGURATION OF TEST SYSTEM	7
4. PRELIMINARY TEST	7
4.1 AC Power line Conducted Emission Test	7
4.2 RADIATED EMISSION TEST	7
5. FINAL RESULT OF MEASURMENT	8
5.1 CONDUCTED EMISSION TEST	8
5.2 RADIATED EMISSION TESTS	10
6. LIST OF TEST EQUIPMENT	11

Report No.: E132R-002



Page 3 of 11 Report No.: E132R-002

Revision History

Issued Report No.	Issued Date	Revisions	Effect Section
E132R-002	February 04, 2013	Initial Issue	All

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)



Page 4 of 11 Report No.: E132R-002

1. VERIFICATION OF COMPLIANCE

APPLICANT : SD Biosensor, Inc.

ADDRESS : C-4th&5th Floor Digital Empire Building 980-3, Yeongtong-dong, Yeongtong-gu, Suwon-si,

Kyonggi-do, Korea

CONTACT PERSON : Kim Jae Young / Instrument Development team manager

TELEPHONE NO : +82-31-300-0422

FCC ID : RPJ01GM51

MODEL NO/NAME : 01GM51

SERIAL NUMBER : N/A

DATE : February 04, 2013

EQUIPMENT CLASS	JBP - Part 15 Class B Computing Device Peripheral			
E.U.T. DESCRIPTION	SD GlucoNavii® Voice N Blood Glucose Meter - Unintentional			
	Radiator			
THIS REPORT CONCERNS	ORIGINAL GRANT			
MEASUREMENT PROCEDURES	ANSI C63.4: 2009			
TYPE OF EQUIPMENT TESTED	PRE-PRODUCTION			
KIND OF EQUIPMENT AUTHORIZATION REQUESTED	Certification			
EQUIPMENT WILL BE OPERATED UNDER FCC RULES PART(S)	FCC PART 15, SECTION 15.101			
MODIFICATIONS ON THE EQUIPMENT TO ACHIEVE COMPLIANCE	No			
FINAL TEST WAS CONDUCTED ON	10 m Semi anechoic chamber			

^{-.} The above equipment was tested by ONETECH Corp. for compliance with the requirement set forth in the FCC Rules and Regulations. This said equipment in the configuration described in this report, shows the maximum emission levels emanating from equipment are within the compliance requirements.

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)



Page 5 of 11 Report No.: E132R-002

2. GENERAL INFORMATION

2.1 Product Description

The SD Biosensor, Inc., Model 01GM51 (referred to as the EUT in this report) is a SD GlucoNavii® Voice N Blood Glucose Meter with NFC function and USB port. This report covers for PC peripheral device only and NFC function will be covered by another test report. Product specification described herein was obtained from product data sheet or user's manual.

CHASSIS TYPE	Plastic
LIST OF EACH OSC. OR	0.707
CRY. FREQ.(FREQ.>=1 MHz)	8 MHz
Electrical Rating	Max 3 Vdc, 500 mA
EXTERNAL CONNECTOR	3.5 Pi 3 pole stereo jack, Customized sensor connector

2.2 Model Differences

-. None

2.3 Related Submittal(s) / Grant(s)

-. Original submittal only



Page 6 of 11 Report No.: E132R-002

2.4 Test System Details

The model numbers for all the equipments that were used in the tested system is:

Model	Manufacturer	FCC ID	Description	Connected to
01GM51	SD Biosensor, Inc.	RPJ01GM51	SD GlucoNavii® Voice N Blood Glucose Meter (EUT)	Notebook PC
PP21L	DELL	DoC	Notebook PC	-
ADP-65HB AD	DELTA Electronics. Inc.	N/A	Notebook PC adaptor	Notebook PC
E176FPb	DELL	DoC	Monitor	Notebook PC
KB38211	HP	DoC	Keyboard	Notebook PC
1113	MICROSOFT	DoC	Mouse	Notebook PC

2.5 Cable Description for the Test System

Cable	Shielded	Ferrite Bead	Metal Shell	Length (m)	Connected to
Stereo jack	Y	Both End	Both End	0.5	USB port of Notebook PC

2.6 Test Methodology

Both conducted and radiated testing was performed according to the procedures in ANSI C63.4: 2009. Radiated testing was performed at a distance of 10 m from EUT to the antenna up to 1 GHz.

2.7 Test Facility

The Electromagnetic compatibility measurement facilities are located on at 301-14, Daessangryung-ri, Chowol-eup, Gwangju-si, Gyeonggi-do, 464-862, Korea. The Onetech Corp. has been accredited as a Conformity Assessment Body (CAB) with designation Number, KR0013.

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)



Page 7 of 11 Report No.: E132R-002

3. SYSTEM TEST CONFIGURATION

3.1 Mode of operation during the test

-. The EUT was connected to the Stereo jack port of notebook PC and transmit the data to notebook PC continuously through the USB port during the test.

3.2 Equipment Modifications

-. None

3.3 Configuration of Test System

Line Conducted Test

: The EUT was connected to LISN. All supporting equipments were connected to another LISN. Preliminary Power line Conducted Emission test was performed by using the procedure in ANSI C63.4: 2009 7.3.3 to determine the worse operating conditions.

Radiated Emission Test

: Preliminary radiated emission test was conducted using the procedure in ANSI C63.4: 2009 8.3.1.1 to determine the worse operating conditions. Final radiated emission test was conducted at 10 m semi anechoic chamber.

4. PRELIMINARY TEST

4.1 AC Power line Conducted Emission Test

During Preliminary Test, the following operating mode was investigated

Operation Mode	The Worse operating condition (Please check one only)
Data communication Mode	X

4.2 Radiated Emission Test

During Preliminary Test, the following operating mode was investigated

Operation Mode	The Worse operating condition (Please check one only)
Data communication Mode	X

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)





5. FINAL RESULT OF MEASURMENT

Preliminary test was done in normal operation mode. And the final measurement was selected for the maximized emission level.

5.1 Conducted Emission Test

Humidity Level : <u>34 % R.H.</u> Temperature: <u>22 °C</u>

Limits apply to : FCC CFR 47, PART 15, SUBPART B, SECTION 15.107(a)

Type of Test : <u>CLASS B</u>

Result : PASSED BY 7.20 dB at 2.18 MHz under Peak detector mode

EUT : SD GlucoNavii® Voice N Blood Glucose Meter Date: February 01, 2013

Detector : CISPR Quasi-Peak (6 dB Bandwidth: 9 kHz)

Frequency		Quasi-Pea	k (dBμV)	Margin
(MHz)	Line	Emission level	Q.P Limits	(dB)
0.20	N	54.00	63.60	9.60
0.21	Н	51.90	63.40	11.50
0.27	N	44.70	61.00	16.30
1.22	N	34.90	56.00	21.10
2.14	N	48.70	56.00	7.30
2.18	Н	48.80	56.00	7.20
Frequency		Average (dBμV)		Margin
(MHz)	Line	Emission level	Limits	(dB)
0.20	N	39.60	53.60	14.00
0.21	Н	39.20	53.40	14.20
2.14	N	37.20	46.00	8.80
2.18	Н	37.80	46.00	8.20

Line Conducted Emissions Tabulated Data

Remark: "H": Hot Line, "N": Neutral Line.

Margin (dB) = Limit - Emission Level

 $Emission\ Level\ (dB\mu V) = Reading\ value\ (dB\mu V) + Insertion\ Loss\ of\ LISN\ (dB) + Cable\ loss\ (dB)$

See next page for an overview sweep performed with quasi-peak and average detector.


Tested by: Sung-Woo, Park / Project Engineer


It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)





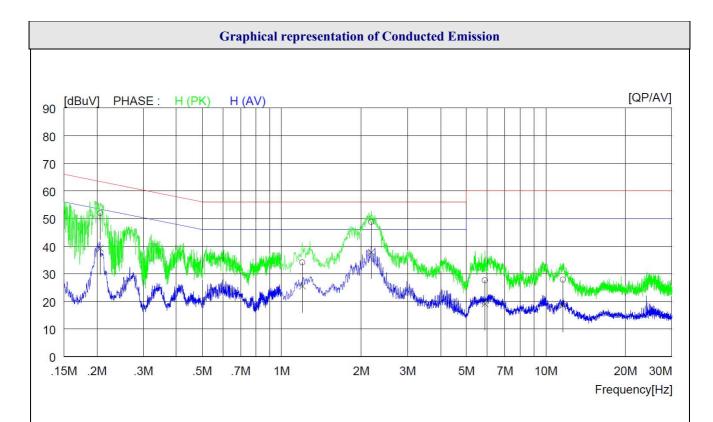

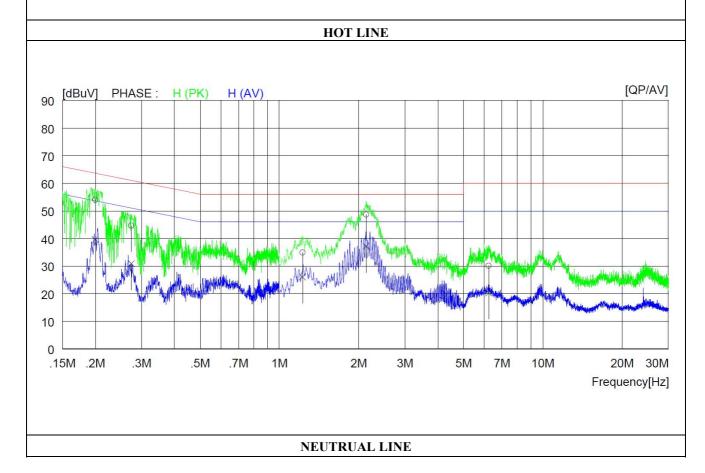

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)

Report No.: E132R-002



Page 10 of 11 Report No.: E132R-002

5.2 Radiated Emission Tests

The following table shows the highest levels of radiated emissions on both polarizations of horizontal and vertical.

Humidity Level : 34 % Temperature: 21 °C

Limits apply to : FCC CFR 47, PART 15, SUBPART B, SECTION 15.109 (g)


Type of Test : CLASS B

Result : PASSED BY 6.20 dB at 47.46 MHz under quasi-peak detector mode

EUT : SD GlucoNavii® Voice N Blood Glucose Meter Date: February 01, 2013

Detector : CISPR Quasi-Peak (6 dB Bandwidth: 120 kHz)

Distance : 10 Meter

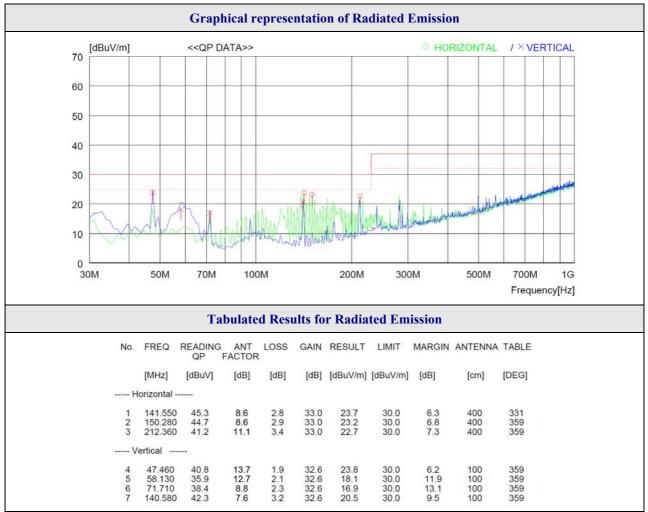

Remark: Margin (dB) = Limit – Result and Result = Reading Peak + Antenna Factor + Loss – Gain Loss and Gain in above table means Cable Loss and Pre-amplifier gain.

Tested by: Sung-Woo, Park / Project Engineer

puf

It should not be reproduced except in full, without the written approval of ONETECH Corp.

EMC-002 (Rev.2)





6. LIST OF TEST EQUIPMENT

No.	Equipment	Manufacturer	Model Name	Serial No.	Last Cal.	Interval Cal.	Used
1.		Rohde & Schwarz	ESCI	101013	Oct. 14, 2012	One Year	
2.	T		ESU	100261	Sep. 24, 2012	One Year	
3.	Test receiver		ESiB26	100296	Apr. 13, 2012	One Year	
4.			ESHS 10	834467/007	Jun. 21, 2012	One Year	
5.		C In -4	310N	312544	May 30, 2012	One Year	
6.	Pre-Amplifier	Sonoma Instrument	310N	312545	May 30, 2012	One Year	
7.		Rohde & Schwarz	SCU 18	10041	Jan. 25, 2012	One Year	
8	TRILOG Broadband		VULB9163	9163-255	Apr. 24, 2012	Two years	
9.	Antenna	Schwarzbeck	VULB9163	9163-420	Mar. 27, 2012	Two years	
10.	Horn Antenna	Schwarzbeck	BBHA9120D	BBHA9120D295	Aug. 23, 2011	Two years	
11.		EMCO	2025/2	9109-1867	May 30, 2012	One Year	
12.	LICN	EMCO	3825/2	9109-1869	May 30, 2012	One Year	
13.	LISN	C -1 11-	NSLK 8126	8126-404	Jun. 11, 2012	One Year	
14.		Schwarzbeck	NSLK 8128	8128-216	Jun. 11, 2012	One Year	
15.	Controller	Innco System	CO2000	619/27030611/L	N/A	N/A	
16.	Turn Table	Innco System	DT3000	930611	N/A	N/A	
17.		T C 4	MA4000-EP	3320611	N/A	N/A	
18.	Antenna Master Innco System		MA4000-EP	3350611	N/A	N/A	
19.	Tripod	EMCO	N/A	N/A	N/A	N/A	

Remark: Mark ■ mean used equipment.

Report No.: E132R-002