


Georgios Kokovidis  
Dräger Medical Systems, Inc.  
6 Tech Drive  
Andover, MA 01810

1-978-379-8148  
george.kokovidis@draeger.com

Ad-hoc mode for TeleSmart M300, patient worn transceiver, can not be implemented by changing a setting from the device user interface. In order to change the operating mode to Ad-hoc mode, one must plug in a special interface programming cable with a custom connector into the UART of the device, and with a laptop, using Hyperterminal, connect to the device. Once connected, you must enter a password to change any of the settings. Ad-hoc is available by the chip set manufacturer, but it is not supported by the firmware resident in the device.

Georgios Kokovidis

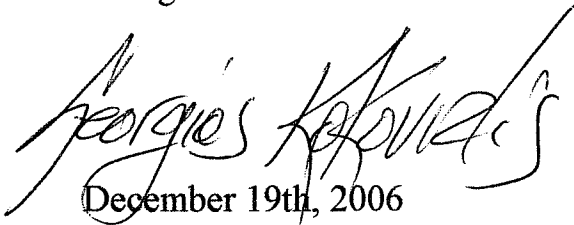

December 19th, 2006

12/19/2006

**Draeger Medical Systems, Inc.**

6 Tech Drive  
Andover, MA 01810 USA  
Phone: (978) 379-8000  
Fax: (978) 379-8332